TRUDELL
HEALTHCARE SOLUTIONS

MEDICAL M&B S.A.S

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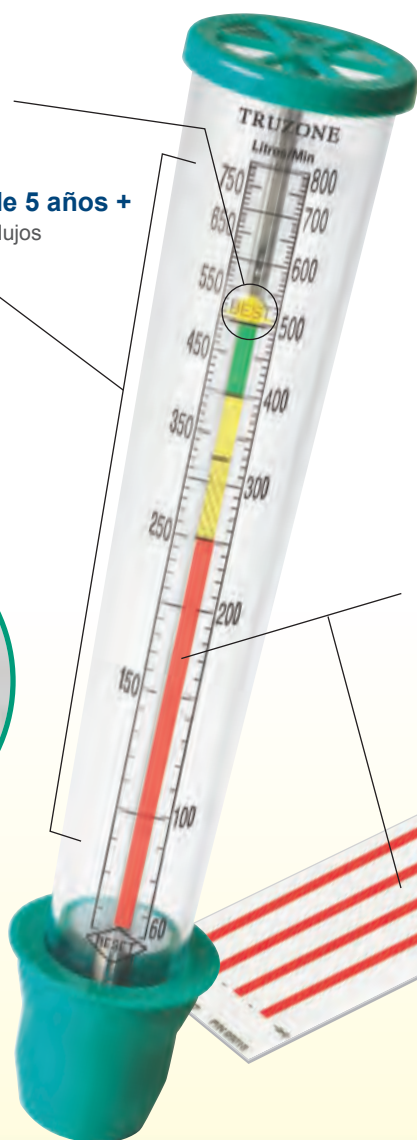
Fácil de usar

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THE LUNG ASSOCIATION

The Lung Association acknowledges that peak flow monitoring is one useful tool to measure airway obstruction.

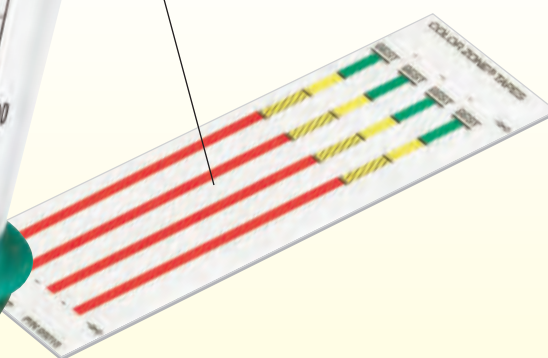
Canadian Asthma Consensus Guidelines recommend the use of a Peak Flow Meter with an Asthma Action Plan for all patient types.¹

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Manufactured by
Monaghan Medical Corporation for:

Trudell Medical International™

725 Third Street, London,
Ontario, Canada, N5V 5G4
Phone: +1-519-455-7060 Fax: +1-519-455-7858
Email: customerservice@trudellmed.com
www.trudellmed.com

www.truzone.ca

TRUZONE® PEAK FLOW METER

1 Gently shake the TRUZONE® Peak Flow Meter (PFM) vertically, until the yellow indicator has been reset, resting within the diamond near the mouthpiece.

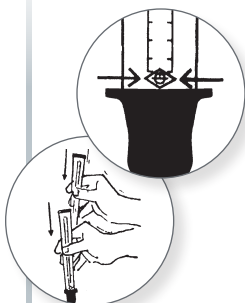
2 Take as deep a breath as possible.

3 Put the TRUZONE® PFM mouthpiece into your mouth, and seal your lips.

4 Blow out as hard and as fast as possible.

5 Look at the number on the scale beside the yellow indicator. Record this number in the space provided in your TRUZONE® PFM Daily Record.

6 Repeat steps 1 through 4 as many times as instructed by your physician or health care professional.

Ask your Physician to do the following:

Use the enclosed ColorZone® Tape to measure your PEF. Place the tape in the middle of the scale so that the arrows below the word "BEST" line up with your best PEF (confirmed by your Physician).

Press the green zone into place and cut the red zone just above the "RESET" diamond. Peel off remaining backing and press into place.

Please refer to the TruZone® PFM package insert for complete instructions and care.

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