Ambu

English

Directions for use

Single Patient Monopolar Needle Electrode (with lead wire)

Indication for use

The Ambu Neuroline Disposable Monopolar needle electrodes are made for muscle activity recording for Electromyography (EMG) applications.

Precautions

Ambu Neuroline Monopolar needle electrodes are for single patient use only. Do not use if individual pouch is damaged or has been opened.

Do not reuse. US: Rx only.

Sterilization/Sterilisation


Sterilized/Sterilised by Ethylene oxide sterilization/sterilisation (ETO).

Packaging

One needle electrode in a pouch.

Directions for use

1.

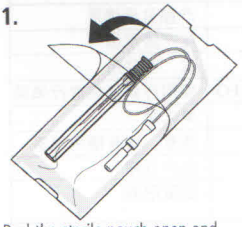Peel the sterile pouch open and remove the single patient needle electrode.

2.

Connect the Ambu Neuroline Monopolar needle electrode with lead wire, that carries a standard 1.5 mm DIN connector, directly to the equipment. Remove the protective tube.

3.

After use, dispose of the needle electrode in a designated container. Do not reuse or re-sterilize/re-sterilise.

Symbol	Indication
REF	Catalog/Catalogue number
LOT	Batch code
	Use by
Sterile EO	Sterilized/sterilised by Ethylene oxide
	Single use
	Attention, see instruction for use
	Date of manufacture
	Temperature limitation
CE 0086	The Monopolar needle electrodes are in conformity with European Medical Directive 93/42/EEC

CE 0086

Ambu® is a registered trademark of Ambu A/S, Denmark.
Ambu A/S is certified according to ISO 9001 and ISO 13485.

Made in Malaysia
Ambu A/S
DK-2750 Ballerup
Tel.: +45 72 25 20 00
Fax: +45 72 25 20 50
www.ambu.com

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Ambu® es una marca registrada de Ambu A/S, Dinamarca.
Ambu A/S dispone de certificado según ISO 9001 e ISO 13485.
Ambu® är ett registrerat varumärke tillhörande Ambu A/S, Danmark.
Ambu A/S är certifierat i enlighet med ISO 9001 och ISO 13485.