
USER MANUAL

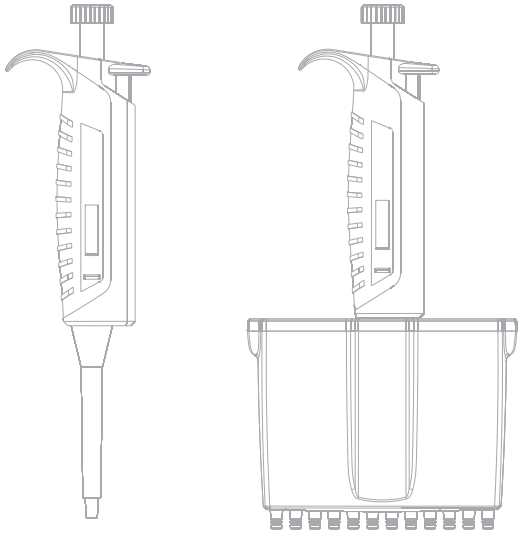SA Series Autoclavable Micropipettes

BOECO Germany

Index

Pipette Description	3-8	Descripción de la pipeta	3-8
Part Description	9	De parte Descripción	9
Pipetting technique	10-12	Técnica de pipeteo	10-12
Calibration and adjustment	13-18	Calibración y ajuste	13-18
Pipetting recommendations	19-20	Recomendaciones de uso	19-20
Maintenance	21	Mantenimiento	21
Cleaning your pipette	22	Limpieza de su pipeta	22
In-house maintenance	23	Mantenimiento interno	23
Autoclaving	24	Autoclavado	24
Recalibration	25	Recalibración	25
Storage	26	Almacenaje	26
Trouble shooting	27-28	Posibles averías	27-28

Pipette Description / Descripción de la pipeta

These PIPETTES are available as fixed or variable, general purpose micropipettes for sampling and dispensing accurate amounts of liquid in micro volumes.

These pipettes operate on the air displacement principle and use disposable tips.

All models are equipped with a built-in Tip Ejector*

They cover the volume range of 0.5 µl to 10000 µl.

* The 10000 µl (10 ml) models BOE 9621111 & BOE 9611111 do have blocked tip ejector. The 10 ml tips have to be removed by hand.

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
Pipetting technique / Técnica de pipeteo

Forward pipetting / Pipeteo directo

1. Press the operating button to the first stop.
2. Dip the tip attached to the pipettes into the solution to a depth of about 1 cm, and slowly release the operating button. Wait for a while then withdraw it from the liquid touching it against the edge of the reservoir to remove excess liquid adhering to outer surface of the tip.
3. Dispense the liquid into the receiving vessel by gently pressing the operating button to the first stop. After a second, press the operating button to the second stop. This will empty the tip completely. Remove the pipette from the vessel sliding it up the wall of the vessel.
4. Release the operating button to the ready position.

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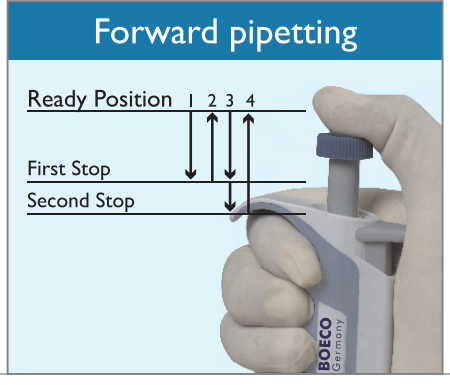
Pipetting technique / Técnica de pipeteo

Reverse pipetting / Pipeteado inverso

1. Press the operating button to the second stop.
2. Dip the tip attached to the pipette into the solution to a depth of about 1 cm, and slowly release the operating button. This action will fill the tip with a volume that is larger than the set volume. Wait 1-2 seconds and withdraw the tip from the liquid touching it against the edge of the reservoir to remove excess liquid.
3. Dispense the liquid into the receiving vessel by pressing the operating button gently and steadily to the first stop. This volume is equal to the set volume. Hold button in this position. Some liquid will remain in the tip, which should not be dispensed.
4. The liquid remaining in the tip can be pipetted back into the original solution or disposed together with the tip.

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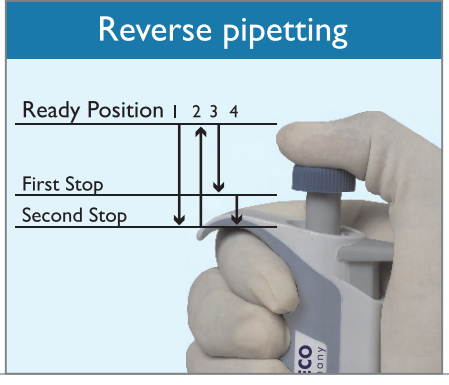Pipetting technique / Técnica de pipeteo

Reverse pipetting / Pipeteado inverso

5. Release the operating button to the ready position.

 Note : Reverse pipetting technique is recommended for viscous solutions, solution having tendency to foam or dispensing very small volumes.

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All pipettes are been quality tested according to ISO8655/DIN12650. The quality control according to ISO8655/DIN 12650 involves gravimetric testing of each pipette with distilled water (quality 3, DIN ISO3696) at 25° C ± 1°C in multi channel pipette, each channel is considered as an individual pipette and checked accordingly.

All pipettes are been calibrated ISO/IEC 17025 accredited laboratory.

Each pipette is calibrated, inspected and validated by qualified technicians according to defined quality system.

Pipettes are constructed to permit re-adjustment for other pipetting techniques or liquids of different temperature and viscosity.

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An analytical balance must be used. The balance selection depends upon selected model of the pipette and sensitivity of balance reading.

Test liquid: Water, distilled or deionized, grade 3 water conforming ISO 3696.

Calibration should be carried out in a draft-free room at a constant (±0.5°C) temperature of water, pipette and air between 15°C to 30°C. The relative humidity must be above 50%, especially with volumes under 50µl, the air humidity should be as high as possible to reduce the effect of evaporation loss.

Special accessories for analytical balance, such as the evaporation trap are recommended for the calibration of volumes under 50µl.

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The pipette is checked at maximum volume, at 50% of maximum volume and at minimum or 10% of maximum volume, whichever is higher.

- A new tip is first pre-wetted 3-5 times and a series of ten pipetting is done at each volume.
- Use of forward pipetting technique is recommended.
- Calculate the inaccuracy and imprecision for all three volumes as per EN ISO 8655/ DIN 12650 standards on the basis of the following calculation.

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Mean Value $\overline{V} = \overline{X} \cdot Z$

Mean Value $\overline{X} = \frac{\sum X_i}{n}$

Xi = Balance Reading, n = number or reading,

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$$A\% = \frac{\bar{V} - Vo}{Vo} 100$$

$\bar{V}$ = Mean Value

Vo = Particular volume at which readings are taken

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$$S = \sqrt{\frac{\sum_{i=1}^{n} (Vi - \overline{V})}{n = 1}}$$

S = Standard Deviation
V = MeanValue
n = number of readings

$$CV\% = \frac{100 \cdot S}{V}$$

- Compare the results to the limits in the subsequent tables

$$S = \sqrt{\frac{\sum_{i=1}^{n} (Vi - \overline{V})}{n = 1}}$$

S = Standard Deviation
V = MeanValue
n = number of readings

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Pipetting recommendations / Recomendaciones de uso

- Aspirate liquid into the pipette only when a tip is attached to its tipcone.
- While pipetting, the pipettor should be vertically straight and tip should be dipped only a few millimeters into the liquid.
- Prerinsing of tip 5 times with the liquid to be dispensed is recommended. This is important especially when dispensing liquids which have a viscosity and density different from water.
- Always control the push button movements with the thumb for consistency.
- Allow liquids, Tips, and pipettes to equilibrate to the ambient temperature.
- Presinse the tip several times before use when pipetting liquids at temperature different from ambient.
- Wipe the tip only if there is liquid on the outside of the tip, being careful to avoid touching tip's orifice.

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Pipetting recommendations / Recomendaciones de uso

- Don't keep pipette in your hand while not working to avoid transferring body heat.
- Use the correct pipette tip designed for use with the particular pipette
- Select the correct pipetting technique (e.g. Reverse, forward etc.) depending on the nature of the liquid.
- Using excessive force to turn the push button outside the range specified for it may jam the mechanism and damage the pipettor.

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Maintenance / Mantenimiento

To maintain the best results from your pipettor, each unit should be checked every day for cleanliness. Particular attention should be paid to the tip cone(s).

This pipettor has been designed for easy in-house service. However, also provides complete repair and calibration service. Please return your pipettor to your local distributor for repair or calibration. Before returning please make sure that it is free from all contamination.

Check the performannce of your pipettor regularly e.g every 3 months and always after in-house service or maintenance.

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Cleaning your pipette / Limpieza de su pipeta

To clean your pipettor use ethanol and a soft cloth or lint-free tissue. It is recommended to clean the tip cone regularly.

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In-house maintenance / Mantenimiento interno

- Hold down the tip ejector.
- Place the tooth of the opening tool between the tip ejector and the tip ejector collar to release the locking mechanism.
- Carefully release the tip ejector and remove the ejector collar.
- Place the wrench end of the opening tool over the tip cone and turn it anti-clockwise.
- After removing the tip cone, wipe the piston, the O-ring and the tip cone with ethanol and a lint-free cloth.

 Note : Models up to 10µl have a fixed O-ring location inside the tip cone, so do not remove or maintain it.
- Before replacing tip cone, it is recommended to grease the piston sligthtly using the silicone grease provided.
- After reassembling, use the pipettor (without liquid) several times to make sure that the grease is spread evenly.
- Check the pipettor calibration.

 Note : Never disassemble the upper part of the pipette. To avoid losing or damaging fragile parts, reassemble the pipette immediately.

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Autoclaving / Autoclavado

Boeco pipette is completely autoclavable at 2 bar pressure 121°c temperature for 20 minutes exposure time according to DIN

- Autoclave the complete pipette without any disassembling.
- After autoclaving, allow pipette to completely cool and dry.

If the pipette is autoclaved frequently, the piston and springs should be greased with the lubricant supplied with each pipette to maintain smooth movement.

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Recalibration / Recalibración

(1) Place the service tool in to the holes at the base of the Push button as shown below. Turn it clockwise to increase & anticlockwise to decrease the volume.

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Storage / Almacenaje

When not in use it is recommended that your pipettor is stored in a vertical position.

Leaving pipette on it side can cause liquids to like in to the body of the pipette and cause corrosion.

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Trouble shooting

Trouble	Possible Reason	Correction
Droplets left inside the tip	Unsuitable tip, Non-Uniform wetting of the plastic	Use new tip
Leakage or Pipetted volime too small	Tip holder (cone) scratched or damaged	Change the tip cone
	Organic Solvent as liquid.	Aspirate & discard the organic Solvent several times before actual pipetting by the same tip.
	Tip incorrectly attached	Attach firmly.
	Unsuitable tip.	Use new tip.
	Foreign particles between tip and tip cone.	Clean the tip cone.
	Insufficient amount of grease on piston and O-ring.	Clean & grease O-ring and piston.
	O-ring not correctly positioned or damaged.	Change the O-ring.

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Trouble shooting

Trouble	Possible Reason	Correction
Inaccuracies	Incorrect operation.	Follow instructions carefully.
	Calibration altered.	Recalibrate according to instruction.
	Unsuitable for the particular liquid pipetting technique.	Use correct pipetting technique.
	Instrument damaged.	Send for repair.
Push button jammed or moves erratically	Piston contamined.	Clean & grease O-ring and piston
	Penetration of solvent vapours.	
Tip Ejector jammed or moves erratically.	Tip cone contaminated from outside.	Remove ejector collar and clean tip cone's outer surface with ethanol.
Volume setting is not properly click stopped.	Click stop mechanism damaged.	Send for repair.
Push Button dose not turn for volume setting.	Use of excessive force beyond the range of pipette	Send for repair.

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