flamenco™

spinal fixation system

Optimale Implantatauswahl
Optimum implant selection

100 Jahre ulrich medical

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
100 years of reliability and German precision

Spinal systems by ulrich medical® stand for quality "Made in Germany."
They are the successful result of systematic development activities and many years of experience in medical technology. For more than 100 years, we have done our best every day for our customers and for patient health with our proprietary and innovative products.

Leading-edge technology and competence from one source

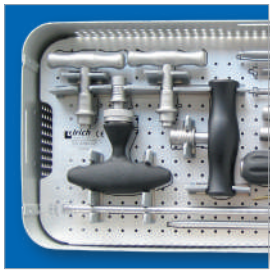		Seite Page	
	Indikationen und Kontraindikationen	4	1
	Indications and contraindications	**4**	
	Sterilisation und Systemkomponenten	5	2
	Sterilisation and System Components	**5**	
	flamenco™ Implantate	6	3
	flamenco™ implants	**6**	
	flamenco™ Instrumente	9	4
	flamenco™ instruments	**9**	
	Operationstechnik	14	5
	Surgical technique	**14**	
	Siebe	22	6
	Trays	**22**	
	Komponenten	23	7
	Components	**23**	

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Intended use

The flamenco™ spinal fixation system is a fixation system for the thoracic and lumbar spine and sacrum. It is composed of pedicle screws, rods, hooks, crosslink assemblies, and connectors for spinal surgery.

Various sizes and styles of the implants are available for optimal adaptation of the pathology of the individual patient. Specific instruments have been designed for use with these implants.

Indications

flamenco™ is an implant system based on pedicle screws indicated for use with skeletally matured patients for rigid fixation and stability who present the acute or chronic instabilities or deformities of the following diagnosis, and receive autogenous bone grafting with above device:

- Degenerative spondylolisthesis with distinctive evidence of neurological impairment
- Vertebral fracture
- Spinal tumor
- Dislocation
- Scoliosis
- Kyphotic deformity
- Discogenic deformity
- Failed previous fusion

Contraindications

- Bone absorption, osteopenia, and/or osteoporosis
- Any active or suspected latent infection of the spine
- Any mental or neuromuscular disorder which might create unacceptable risk of fixation failure or complications post-operatively
- Bone stock abnormalities, or deficiency which can not provide adequate support and/or fixation to the implants
- Pathological obesity
- Open wounds
- Metal sensitivity, documented or suspected
- Pregnancy
- Excessive local inflammation reaction
- Other medical or surgical symptoms that may preclude the potential benefit of spinal implant surgery, such as the presence of tumors, congenital abnormalities, elevation of sedimentation rate unexplainable by other diseases, elevation of white blood count (WBC), or marked left shift in the WBC differential count.

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Sterilization

ulrich medical® suggests the devices to be sterilized via steam sterilization. Please observe the instructions for use for more detailed information.

System Components

Components of the flamenco™ system afford a high extent of versatility and simplicity in the course of spine surgery.

- All key components including screw series, rod series, hook series, crosslink series, are available in a variety of specifications and dimensions for the needs of comprehensive spinal surgery requirements.
- Periphery components such as locking screws, connector blocks are available in complete range to accommodate key components for easy assembly to the users.

Please note the instructions for use for flamenco™ implants and instructions for reprocessing treatment for instruments.

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The small pictures at the bottom of the page show the step-by-step application of the instruments that are used as per the surgical steps on the double page. Pictures with instruments that had been used before are blue-colored.

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Ø 7,5 mm blau, Länge 35, 40, 45, 50, 55 mm

flamenco™ screws, monoaxial

Four diameters in different lengths:

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Ø 8,0 mm bronze, Länge 35, 40, 45, 50, 55 mm

flamenco™ screws, polyaxial

Five diameters in different lengths:

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Ø 6,5 mm grün, Länge 35, 40, 45, 50, 55 mm

Ø 7,5 mm blau, Länge 35, 40, 45, 50, 55 mm

flamenco™ screws, long-arm, monoaxial

Three diameters in different lengths:

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flamenco™ locking screw, silver

Caution: Do not use with screw width Ø 8.0 mm

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flamenco™ rod Ø 5.5 mm
length 40 mm to 350 mm

Additional rods and hybridrods on request.

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flamenco™ rod-rod connector

parallel, Ø 5.5 to 5.5 mm
parallel, Ø 5.5 to 6.0 mm
axial, Ø 5.5 to 5.5 mm
axial, Ø 5.5 to 6.0 mm

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
Stabhaltezange, einstellbar — Rod clamp, adjustable

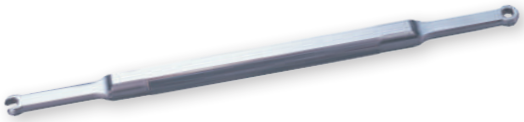CS 8040

Stabrotationsinstrument, SW 5,0 mm — Rotator bar, hex 5.0 mm

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Screwdriver for locking screw, hex 4.0 mm

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***NOTE:** According to the Medical Device Directive 93/42/EEC the torque value of the flamenco™ instrument CS 8049 Torque wrench with T-handle has to be calibrated annually. The nominal torque value is 12.4 Nm ±10%.

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
Approach

With the patient prone on table, the posterior elements are exposed by careful soft tissue dissection.

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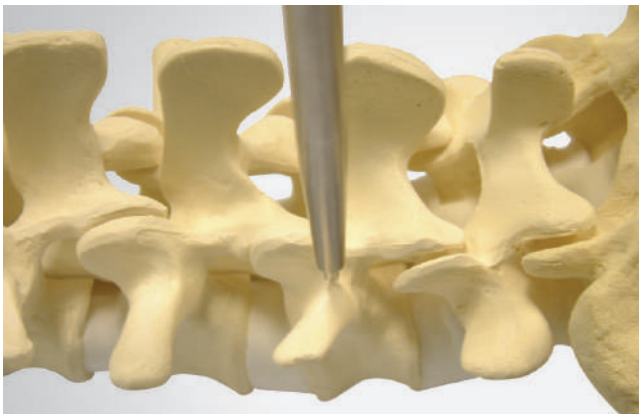
Opening of the Cortical Bone

The cortical bone at the pedicle entry point is opened with the straight awl (CS 8022).

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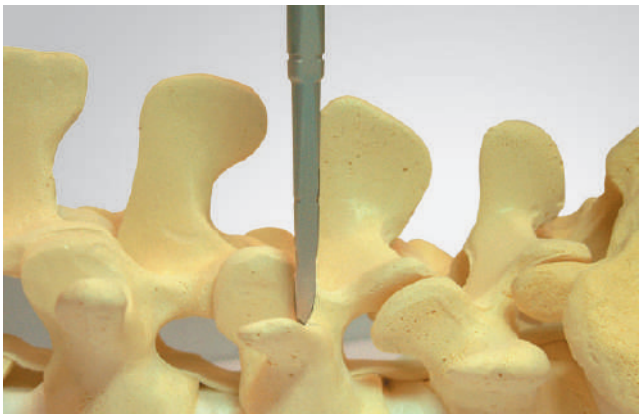
Pedicle Preparation

Using the curved awl (CS 8023), the screw hole is prepared. The curved awl can also be used for determination of the screw length.

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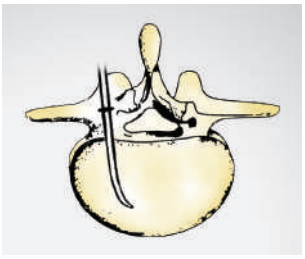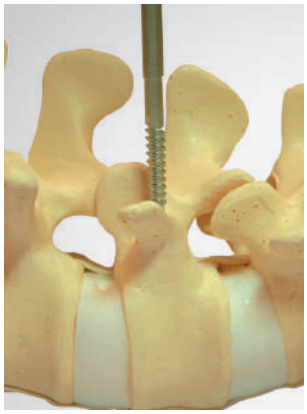
Pedicle Tapping

All flamenco™ screws are self-tapping. If desired, the tap (CS 8026-05, -06, or -07) in combination with the T-handle (CS 8020) can be used to prepare the screw thread in the pedicle hole. The tapped pedicles can be checked again with the ball tip probe.

In the flamenco™ system taps with different diameters are included.

Pedikelschraube **Pedicle screw**	**Instrument** **Instrument**
Ø 5,5 mm Ø 5.5 mm	Gewindeschneider, Ø 5,0 mm (CS 8026-05) Tap, Ø 5.0 mm (CS 8026-05)
Ø 6,5 mm und Ø 7,0 mm Ø 6.5 mm and Ø 7.0 mm	Gewindeschneider, Ø 6,0 mm (CS 8026-06) Tap, 6.0 mm (CS 8026-06)
Ø 7,5 mm und Ø 8,0 mm Ø 7.5 mm and Ø 8.0 mm	Gewindeschneider, Ø 7,0 mm (CS 8026-07) Tap, 7.0 mm (CS 8026-07)

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
Control of the Screw Holes

For lateral C-arm control of the prepared screw holes two different contoured probes (CS 8042 und CS 8043) are available.

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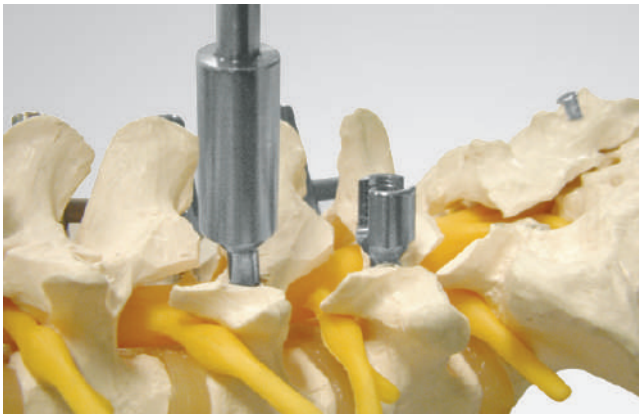
Inserting the Screw

The appropriate screw is mounted and fastened onto the screwdriver.

For protecting the tissue the protection sleeve is used.

Instrument **Instrument**	**Pedikelschraubensystem** **Pedicle Screw System**
Schraubendreher für monoaxiale Schrauben (CS 8028) Screwdriver for monoaxial screws (CS 8028)	Monoaxiale Schrauben und Langkopfschrauben Monoaxial screws and long-arm screws
Schraubendreher für polyaxiale Schrauben (CS 8029) Screwdriver for polyaxial screws (CS 8029)	Polyaxiale Schrauben Polyaxial screws
Schraubendreher für polyaxiale Schrauben Ø 8,0 mm (CS 8030) Screwdriver for polyaxial screws Ø 8.0 mm (CS 8030)	Polyaxiale Schrauben Polyaxial screws
Schraubendreher für polyaxiale Langkopfschrauben (CS 8044-01) Screwdriver for polyaxial long-arm screws (CS 8044-01)	Polyaxiale Langkopfschrauben Polyaxial long-arm screws
Schraubendreher für polyaxiale Langkopfschrauben, für Schrauben Ø 8,0 mm (CS 8044-02) Screwdriver for polyaxial long-arm screws, for screws Ø 8.0 mm (CS 8044-02)	Polyaxiale Langkopfschrauben Polyaxial long-arm screws

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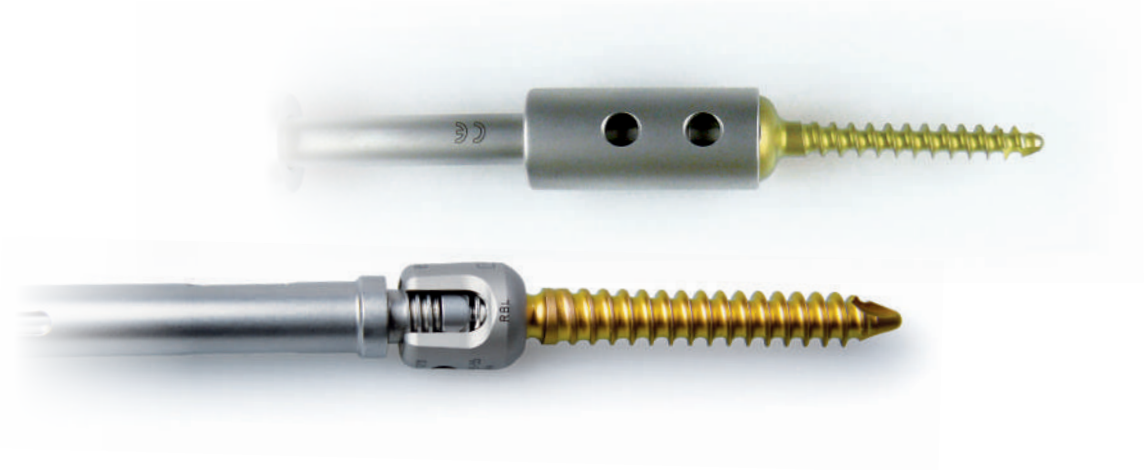The tip of the screw is placed into the pedicle entry point and the screw is inserted using fingertip pressure on the T-handle only.

Great care is taken to maintain screw orientation identified during previous procedure.
When fully seated, the pedicle screws are positioned at depth of 50 % to 80 % of the vertebral body, paralleling the superior end plates. The screwdriver, hex 3.5 mm (CS 8031) is used to adjust the polyaxial screws.

Note: Touch only the rotating sleeve of screwdriver (CS 8029) during screw insertion to prevent loosing of the screwhead fixation.

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Hook Insertion

The hooks are attached to the strategic vertebrae, preoperatively chosen by the surgeon. They can be used isolated or associated to make a clamp. The appropriate inserter (CS 8060 or CS 8061) is used to prepare the hook site.

The desired hook is held with the hook driver (CS 8063) or the hook clamp (CS 8064) and is implanted in the prepared site. If necessary, the hook impactor (CS 8062) can assist the hook insertion.

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
Bending the Rod

The length and the contour of the rod can be determined by the phantom rod (CS 8032-01 or -02).

After choosing the proper length the rods are bent to the appropriate contour by the rod bending forceps (CS 8033). The instrument offers 3 different bending radii. Alternatively, the rod bending forceps, in situ (CS 8034-01 and -02) can be applied for bending.

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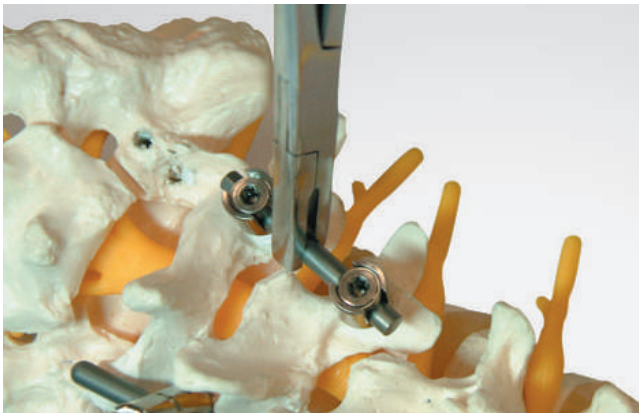Placing the Rod

The rod holder (CS 8035) is used to hold the rod for advancing it into the screw head of the flamenco™ screws.
The rod can easily be advanced in the screw head by using the rod inserter (CS 8036). Place the instrument over the screw head and close the forceps. Via the instruments sleeve the locking screw can safely be fixed in the screw head. Additionally, the rod pusher (CS 8037) or the Rocker (CS 8058) can assist the rod placement.

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
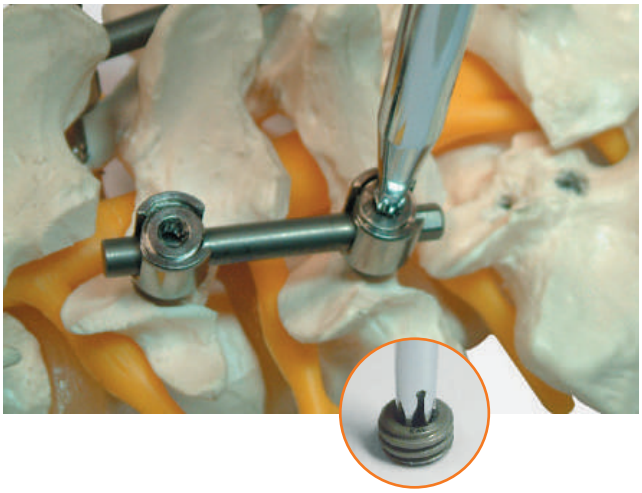
Locking Screw Insertion

Once the rod is seated in all screw heads, the flamenco™ locking screw (CS 8010/CS 8011) in combination with the holder for locking screw (CS 8047) is pre-fixed.
Caution: Do not use flamenco™ locking screw CS 8010 with screw width Ø 8.0 mm!

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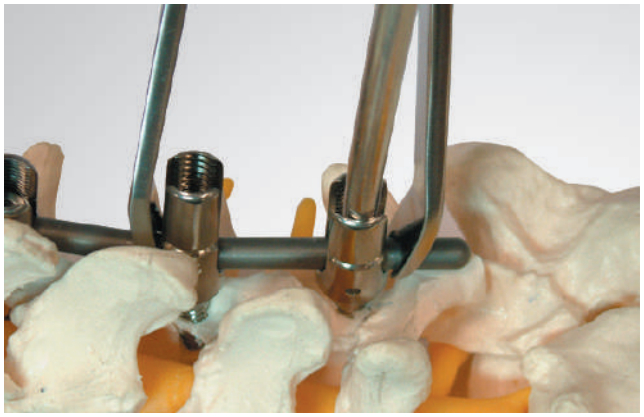
Distraction and Compression

Subject to the clinical situation, distraction and compression are carried out at both sides at the same time by using the distraction forceps, angled (CS 8051) and the compression forceps, angled (CS 8050).

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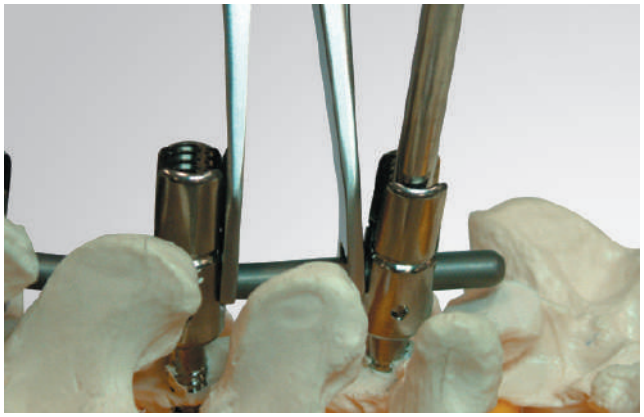Final Tightening

The locking screws at the cranial and caudal ends of the instrumentation are subsequently tightened on both sides, followed by the other adjacent locking screws.

Once compression/distraction and any other disc preparation is completed, the screwdriver for locking screws, hex 4.0 mm (CS 8046) and the torque wrench with T-handle (CS 8049) are assembled. The torque wrench is positioned on the locking screw through the cannulated body of the counter part for screwdriver CS 8046 (CS 8053). Each locking screw is then final tightened sequentially with care.

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Cutting off

The long arm cutting forceps (CS 8045) is used to remove the remaining tabs of the long arm screws.

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Placing the Crosslink

If desired, a crosslink (CS 8013-04-xxx) can be attached to secure the construct stability against rotation and lateral bending.

After selecting the appropriate length of the crosslink, the component is positioned on the instrumentation. The locking screw for crosslink is tightened with the screwdriver for crosslink (CS 8054).

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Implant Removal

Removal of the flamenco Spinal System is performed by reversing the order of the implant procedure. The screwdriver for locking screws, hex 4.0 mm (CS 8046) must first be used to remove the locking screws. The T-handle (CS 8020) in combination with the screwdriver for monoaxial screws (CS 8028) must be used to remove the monoaxial screws.


The T-handle in combination with the screwdriver for polyaxial screws (CS 8029 or CS 8030) must be used to remove the polyaxial screws. The crosslink implants can be removed with the screwdriver for crosslinks (CS 8054).

For emergencies we offer you the explantation system uRevision Set.

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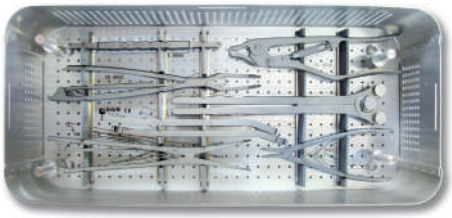
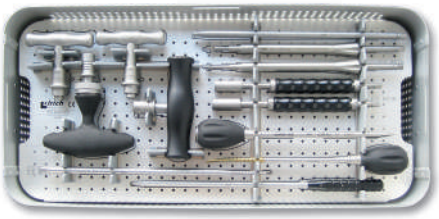
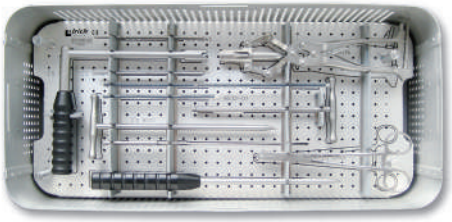
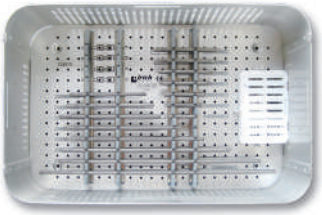
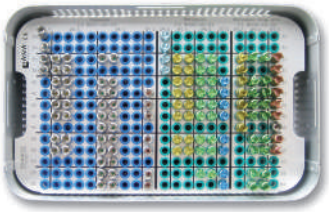
flamenco™ tray for hooks
and instruments

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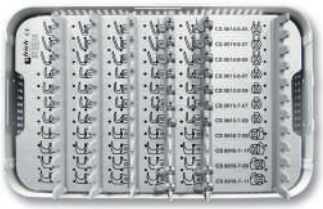
flamenco™ & cosmicMIA™

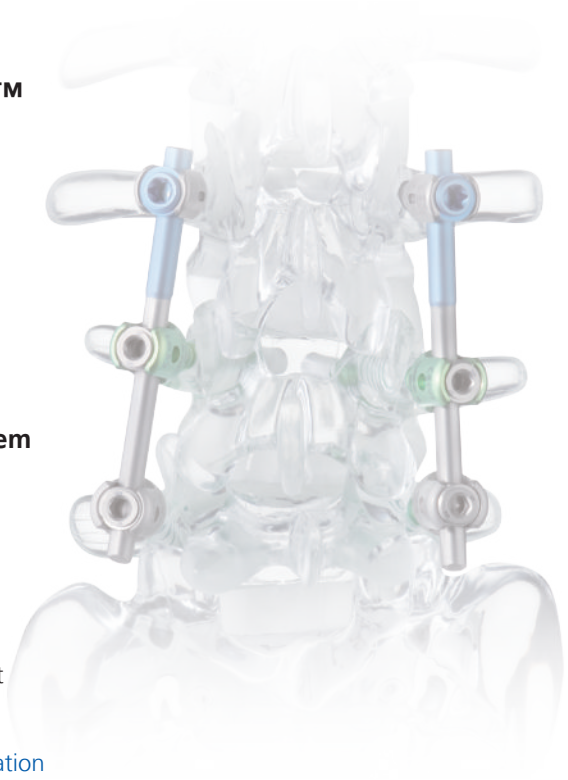cosmicMIA™

dynamic minimally invasive system

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- Posterior system for dynamic stabilization of the thoracic and lumbar spine
- Dynamic screw
- Non-fusion technique
- Preservation of mobility in the segment
- Minimally invasive instrumentation

CE 0123

WS 8001 R6/2014-01

ulrich GmbH & Co. KG I Buchbrunnenweg 12 I 89081 Ulm I Germany
Telefon/Phone: +49 (0)731 9654-0 I Fax national/international: +49 (0)731 9654-2705/2805
spine@ulrichmedical.com I www.ulrichmedical.com

ulrich medical