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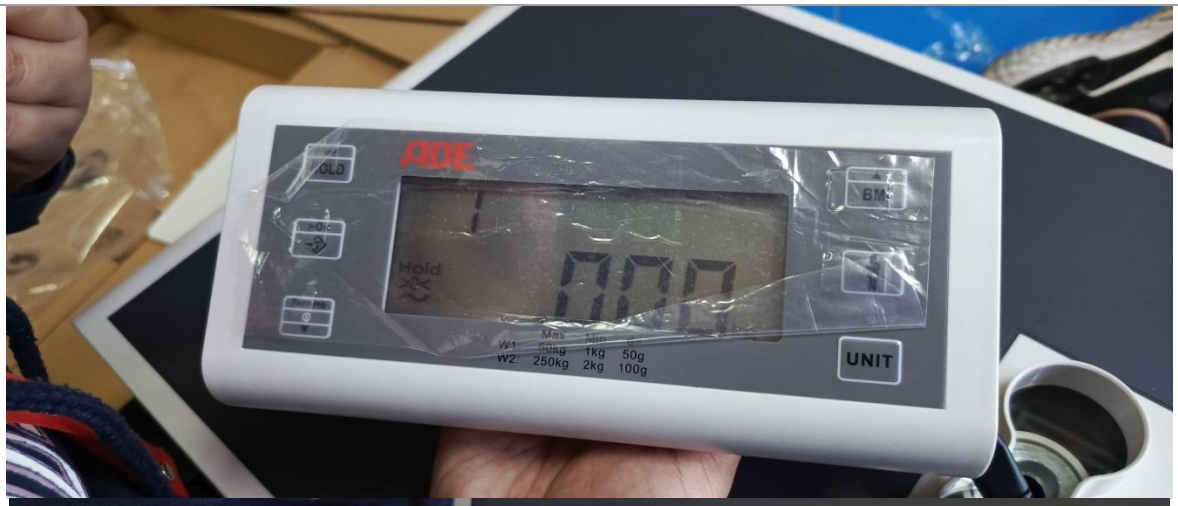
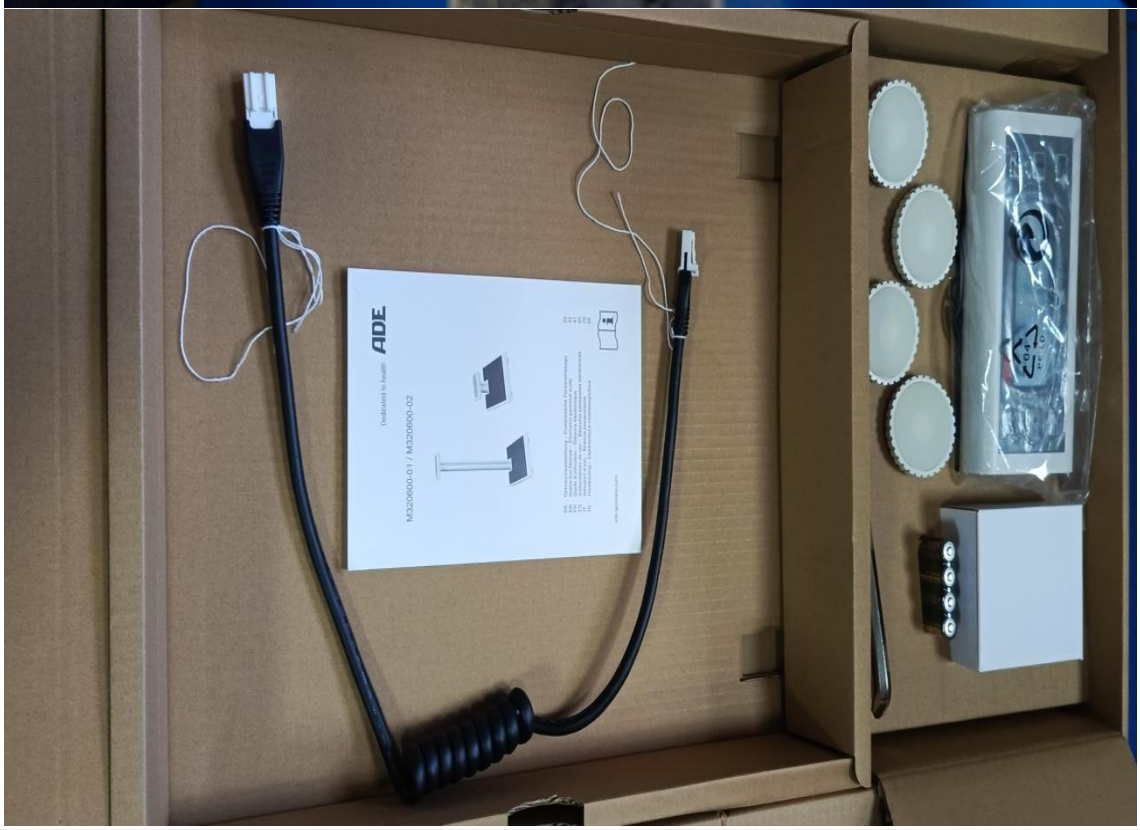
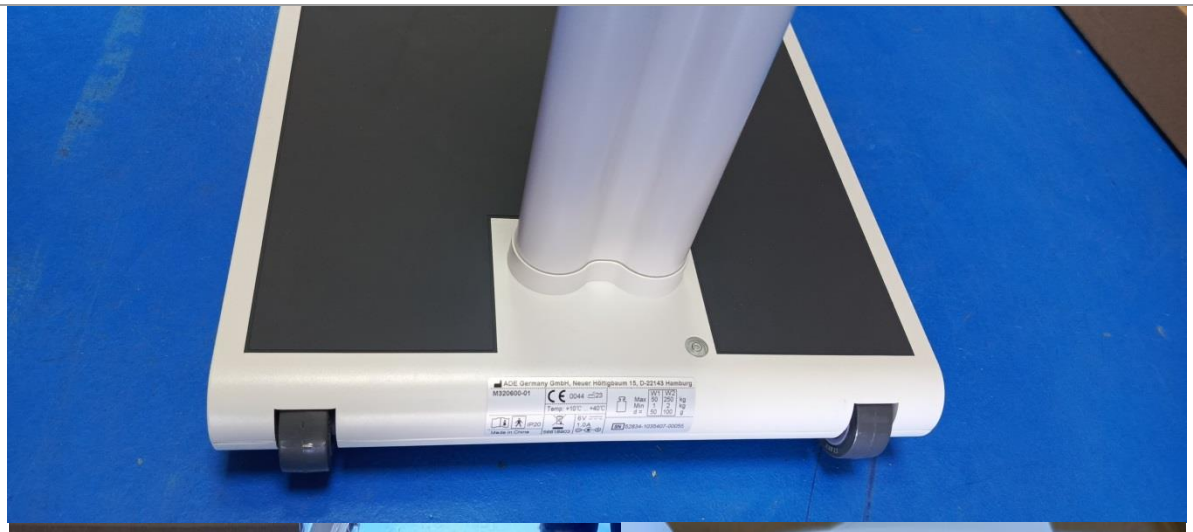
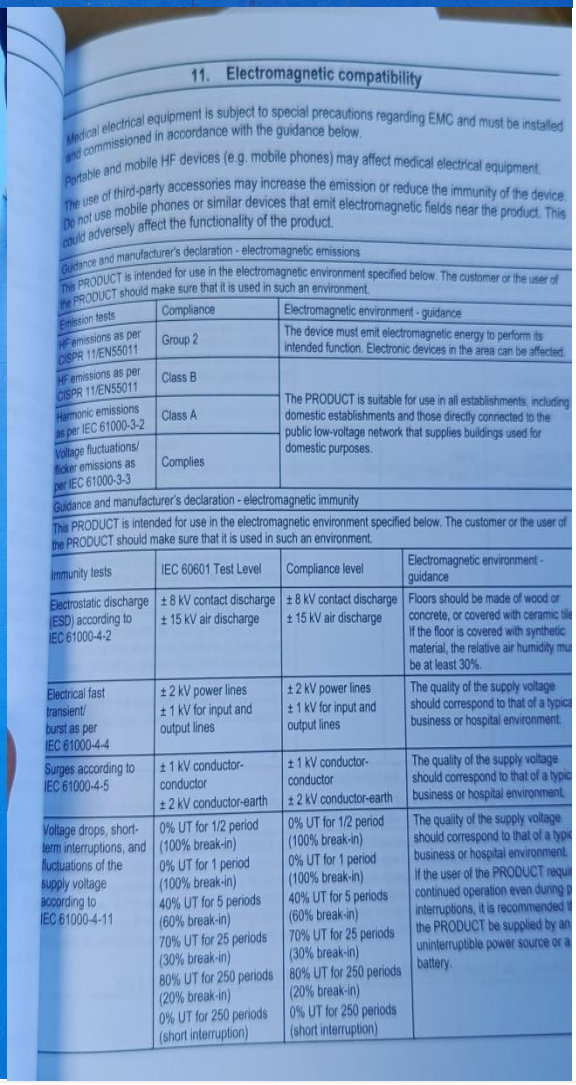
ADE

11. Electromagnetic compatibility

Medical electrical equipment is subject to special precautions regarding EMC and must be installed and commissioned in accordance with the guidance below.

Portable and mobile HF devices (e.g. mobile phones) may affect medical electrical equipment.

The use of third-party accessories may increase the emission or reduce the immunity of the device.

Do not use mobile phones or similar devices that emit electromagnetic fields near the product. This could adversely affect the functionality of the product.

Guidance and manufacturer's declaration - electromagnetic emissions

This PRODUCT is intended for use in the electromagnetic environment specified below. The customer or the user of the PRODUCT should make sure that it is used in such an environment.

Emission tests	Compliance	Electromagnetic environment - guidance
HF emissions as per CISPR 11/EN55011	Group 2	The device must emit electromagnetic energy to perform its intended function. Electronic devices in the area can be affected.
HF emissions as per CISPR 11/EN55011	Class B	The PRODUCT is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage network that supplies buildings used for domestic purposes.
Harmonic emissions as per IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions as per IEC 61000-3-3	Complies	

Guidance and manufacturer's declaration - electromagnetic immunity

This PRODUCT is intended for use in the electromagnetic environment specified below. The customer or the user of the PRODUCT should make sure that it is used in such an environment.

Immunity tests	IEC 60601 Test Level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) according to IEC 61000-4-2	± 8 kV contact discharge ± 15 kV air discharge	± 8 kV contact discharge ± 15 kV air discharge	Floors should be made of wood or concrete, or covered with ceramic tile. If the floor is covered with synthetic material, the relative air humidity must be at least 30%.
Electrical fast transient/ burst as per IEC 61000-4-4	± 2 kV power lines ± 1 kV for input and output lines	± 2 kV power lines ± 1 kV for input and output lines	The quality of the supply voltage should correspond to that of a typical business or hospital environment.
Surges according to IEC 61000-4-5	± 1 kV conductor-conductor ± 2 kV conductor-earth	± 1 kV conductor-conductor ± 2 kV conductor-earth	The quality of the supply voltage should correspond to that of a typical business or hospital environment.
Voltage drops, short-term interruptions, and fluctuations of the supply voltage according to IEC 61000-4-11	0% UT for 1/2 period (100% break-in) 0% UT for 1 period (100% break-in) 40% UT for 5 periods (60% break-in) 70% UT for 25 periods (30% break-in) 80% UT for 250 periods (20% break-in) 0% UT for 250 periods (short interruption)	0% UT for 1/2 period (100% break-in) 0% UT for 1 period (100% break-in) 40% UT for 5 periods (60% break-in) 70% UT for 25 periods (30% break-in) 80% UT for 250 periods (20% break-in) 0% UT for 250 periods (short interruption)	The quality of the supply voltage should correspond to that of a typical business or hospital environment. If the user of the PRODUCT requires continued operation even during power interruptions, it is recommended that the PRODUCT be supplied by an uninterruptible power source or a battery.

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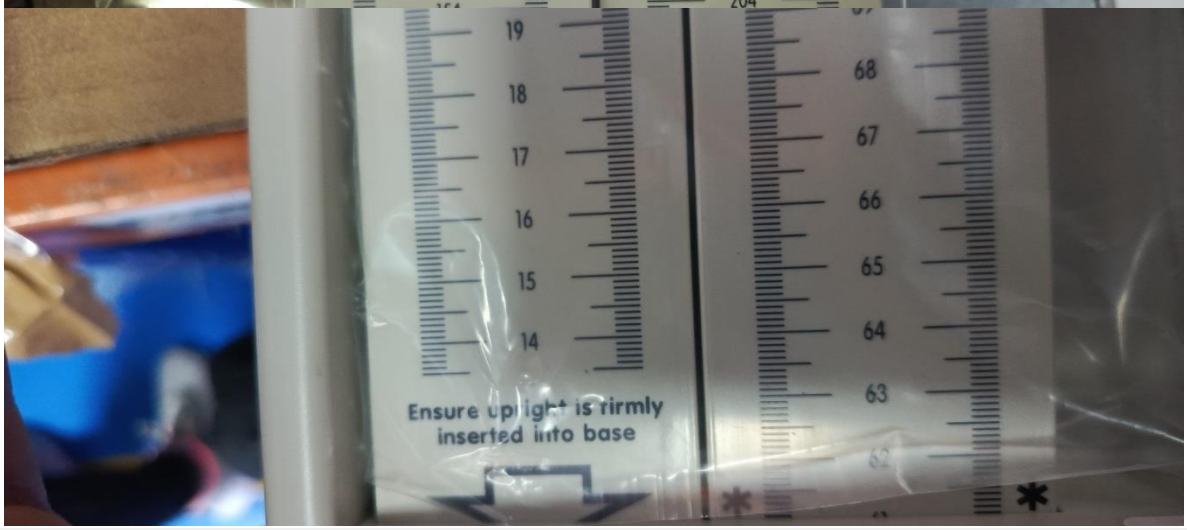
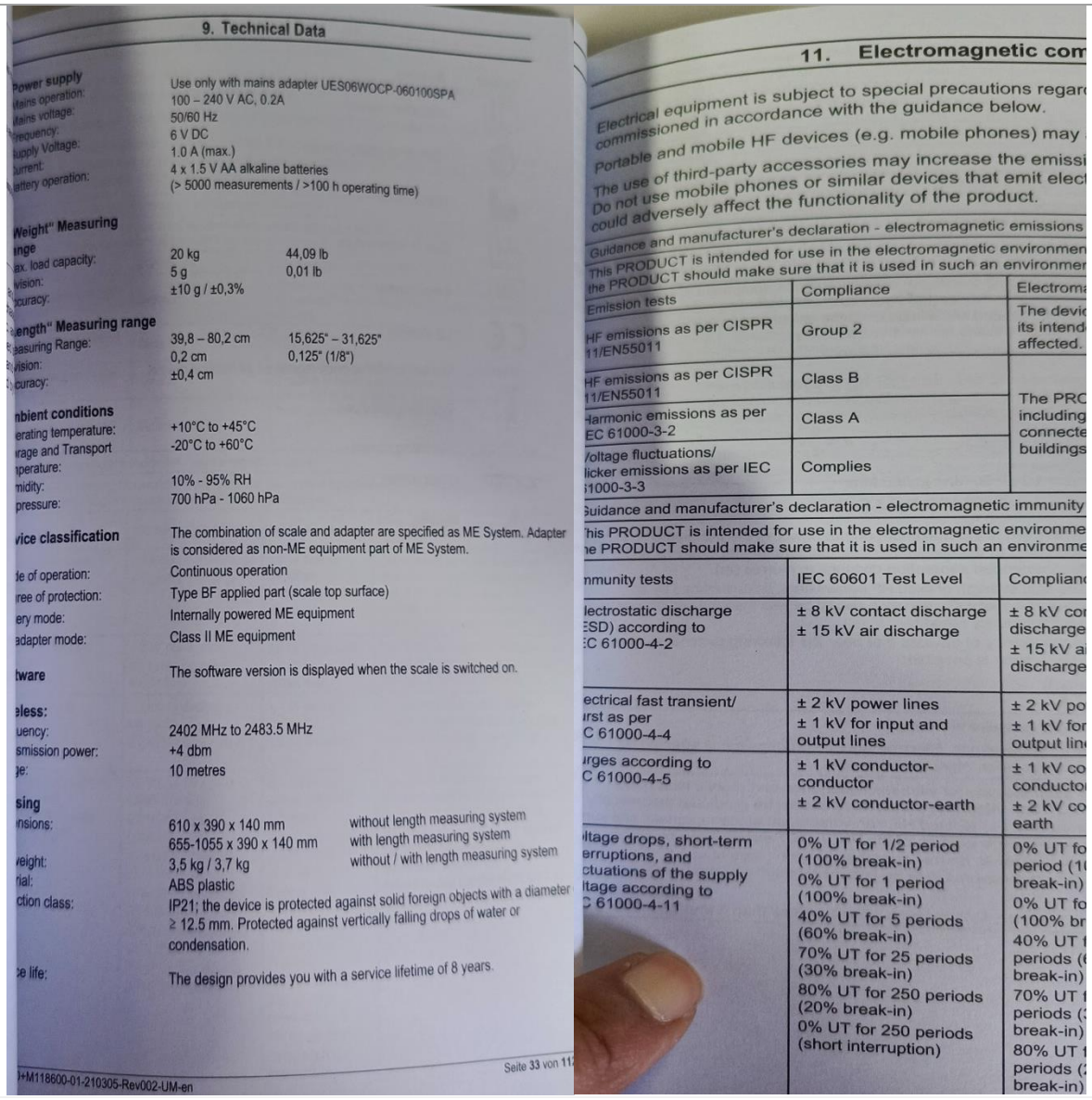
9. Technical Data

Power supply		
Mains operation:	Use only with mains adapter UES06WOCP-060100SPA	
Mains voltage:	100 – 240 V AC, 0.2A	
Frequency:	50/60 Hz	
Supply Voltage:	6 V DC	
Current:	1.0 A (max.)	
Battery operation:	4 x 1.5 V AA alkaline batteries (> 5000 measurements / >100 h operating time)	
Weight" Measuring		
nge		
ax. load capacity:	20 kg	44,09 lb
vision:	5 g	0,01 lb
ccuracy:	±10 g / ±0,3%	
ength" Measuring range		
easuring Range:	39,8 – 80,2 cm	15,625" – 31,625"
vision:	0,2 cm	0,125" (1/8")
curacy:	±0,4 cm	
nbient conditions		
erating temperature:	+10°C to +45°C	
rage and Transport perature:	-20°C to +60°C	
midity:	10% - 95% RH	
pressure:	700 hPa - 1060 hPa	
vice classification	The combination of scale and adapter are specified as ME System. Adapter is considered as non-ME equipment part of ME System.	
e of operation:	Continuous operation	
ree of protection:	Type BF applied part (scale top surface)	
ery mode:	Internally powered ME equipment	
adapter mode:	Class II ME equipment	
tware	The software version is displayed when the scale is switched on.	
eless:		
uency:	2402 MHz to 2483.5 MHz	
smission power:	+4 dbm	
e:	10 metres	
sing		
nsions:	610 x 390 x 140 mm	without length measuring system
	655-1055 x 390 x 140 mm	with length measuring system
eight:	3,5 kg / 3,7 kg	without / with length measuring system
rial:	ABS plastic	
ction class:	IP21; the device is protected against solid foreign objects with a diameter ≥ 12.5 mm. Protected against vertically falling drops of water or condensation.	
e life:	The design provides you with a service lifetime of 8 years.	

11. Electromagnetic com

Electrical equipment is subject to special precautions regar commissioned in accordance with the guidance below.

Portable and mobile HF devices (e.g. mobile phones) may

The use of third-party accessories may increase the emissi

Do not use mobile phones or similar devices that emit elec could adversely affect the functionality of the product.

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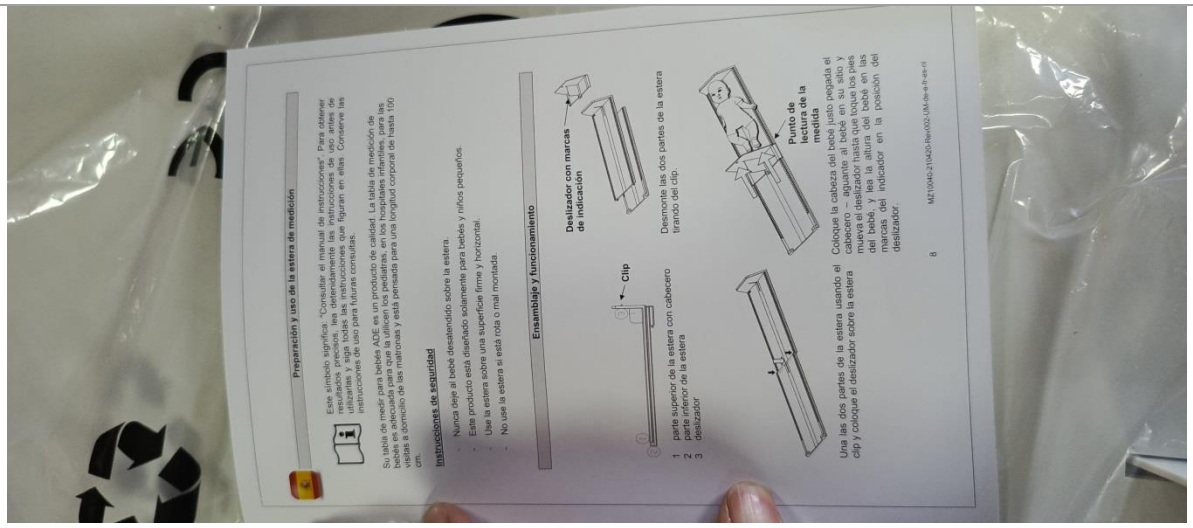
This PRODUCT is intended for use in the electromagnetic environmen the PRODUCT should make sure that it is used in such an environmer

Emission tests	Compliance	Electroma
HF emissions as per CISPR 11/EN55011	Group 2	The devi its intend affected.
HF emissions as per CISPR 11/EN55011	Class B	The PRC including connecte buildings
Harmonic emissions as per IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions as per IEC 61000-3-3	Complies	

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This PRODUCT is intended for use in the electromagnetic environme he PRODUCT should make sure that it is used in such an environme

nmunity tests	IEC 60601 Test Level	Complian
lectrostatic discharge ESD) according to C 61000-4-2	± 8 kV contact discharge ± 15 kV air discharge	± 8 kV co discharge ± 15 kV a discharge
ectrical fast transient/ rst as per C 61000-4-4	± 2 kV power lines ± 1 kV for input and output lines	± 2 kV po ± 1 kV for output lin
rges according to C 61000-4-5	± 1 kV conductor-conductor ± 2 kV conductor-earth	± 1 kV co conducto ± 2 kV co earth
ltage drops, short-term erruptions, and ctuations of the supply ltage according to C 61000-4-11	0% UT for 1/2 period (100% break-in) 0% UT for 1 period (100% break-in) 40% UT for 5 periods (60% break-in) 70% UT for 25 periods (30% break-in) 80% UT for 250 periods (20% break-in) 0% UT for 250 periods (short interruption)	0% UT fo period (1 break-in) 0% UT fo (100% br 40% UT periods (break-in) 70% UT periods (break-in) 80% UT periods (break-in)

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