FORMATO
FICHAS TÉCNICA

Bogotá D.C, 05/06/2026

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Ultrasonic Cleaner / UA22MFDN

Technology

The cavitation phenomenon consists of formation and growth of millions of microscopic bubbles under a low-pressure stage, produced into the ultrasonic bath due to a high-power transducer 37kHz. These bubbles break in all directions and attack every surface in contact with water, erasing dirt, impurities and pollutants when contacting with the materials or pieces into the bath. This phenomenon allows cleaning all kind of piece without need of disassembling, due to it is taking part in all surface in contact with the liquid. The heat allows accelerate this clearing process.

Avant Series

- Entails all advantages of the Eco series
- LCD touch screen for simple and intuitive operation
- Features a digital timer with up to 199 min
- Working power adjustable between 40% and 100%
- Heating ability between 20℃ and 80℃
- Comes with a tray and lid
- Stainless steel 304 tank and paint coated case, meet the requirements of antiseptic and hygienic lab environment

- Independent degas function
- Automatically save last parameters and recalled when unit is turned back on
- Non-slip feet and flexible surface contact for lower ultrasonic vibration
- Safety handles make unit portable
- Maintain ultrasonic efficiency at different temperatures
- Frequency overload and thermal protection

Application Areas

- Laboratories: laboratory material, precision instruments, pipettes, sieves, etc.
- Medicine: instrumentation in general, surgical material, etc.
- Odontology: dental prosthesis, instrumentation, etc.
- Optics: instrumentation, optical, frames, etc.
- Industry: printed circuits, Electronic components, etc.
- Jewellery: watches, jewels, etc.
- Automation: Injector clearing.
- Degasification or liquid dissolution.
- Compact substances degermation

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WIGGENS China
Room 303, Hall C, Office Building M8, No.1 Jiuxianqiao East Road, Chaoyang District, Beijing 100015, China
Tel: +86 400-809-2068
Fax: +86 400-809-2068-112
Email: info@ wiggens.com service@wiggens.com
Website: www.wiggens.com


Multi-Tube Vortex Mixer

Features:

* Mixing mode changes quickly.
* LCD Integrated display.
* High speed accuracy.
* Suitable for a variety of test tube racks.
* The appearance and operation interface are simple and beautiful.

Picture:

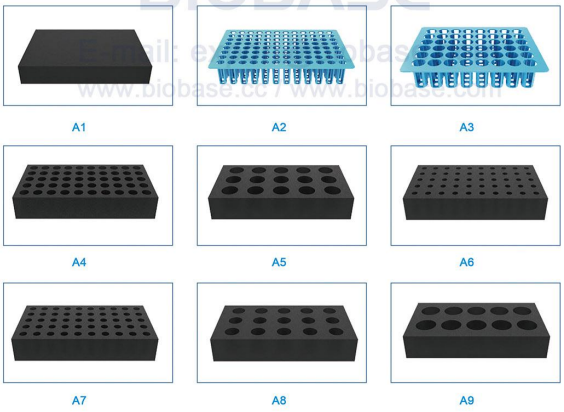A1 A2 A3

A4 A5 A6

A7 A8 A9

BIOBASE
E-mail: export@biobase.com
www.biobase.cc / www.biobase.com

Technical Parameters:

Model	BK-VX2800
Speed Range	300~2800rpm
Speed Regulation Accuracy	±1%
Amplitude	3mm
Time Setting Range	1s~99h59min
Pulse Interval Timing Range	1s~99s
Pulse Operation Timing Range	1s~999s
Max Loading	5.5kg
Consumption	70W
Power Supply	AC100~240V, 50/60Hz
Top Cover Size	310*210mm
External Size (L*W*H)	420*305*465mm
Net Weight	20kg
Package Size (W*D*H)	492*378*520mm
Gross Weight	23kg

Accessories:

Model	Specification	Compatible Tube Type	Hole NO.	Size
Standard				
A1	Replaceable tray mat	/	/	305*204*50mm
A2	Φ17mm plastic test tube rack	96*10ml virus sampling tube can be placed	96	294*200*50mm
A3	Φ26mm plastic test tube rack	Can put 20 mixed 1 virus sampling tube rack	48	280*214*53mm
A4	Φ16mm plastic test tube rack	Can put 10ml/15ml centrifugal tube	50	245*132*45mm
A5	Φ29mm plastic test tube rack	50ml centrifuge tube	15	245*132*45mm
Optional				
A6	Φ10mm(Φ12mm) sponge test tube rack	/	50	245*132*45mm
A7	Φ13mm sponge test tube rack	/	50	245*132*45mm
A8	Φ13mm sponge test tube rack	/	15	245*132*45mm
A9	Φ13mm sponge test tube rack	/	10	245*132*45mm

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Order NO.: W6060050

021

Specific parameters will be adjusted according to product upgrades, without notice.

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Evaporate laboratory solvents in maximum safety

Is this information helpful? 0

Copyright © 2008 WIGGENS CO., LTD. All rights reserved.

Chemical Resistant Diaphragm Pumps C300

Application Areas

> Chemical and petrochemical Industry
> Pharma Industry
> Filtration processes
> Vacuum distillation
> Rotary evaporation
> Vacuum and centrifugal concentration
> Solid phase extraction
> Conventional drying and gel drying
> Advanced substitute for water-jet pumps

C series' construction with PTFE makes them very resistant to chemical vapors from inlet to exhaust and very tolerant to condensates. Pumping chambers are hermetically ensuring long lifetimes of mechanical parts.

Most importantly, diaphragm pumps are oil-free, for vastly reduced service demands compared with oil sealed pumps. They eliminate the cost of water and its contamination well-known from water-jet aspirators, and the waste-oil disposal of rotary vane pumps.
They are chemical resistant diaphragm vacuum pumps which can be widely used to hard acidic, basic and solvent vapors by utilizing corrosion proof PTFE on all the wetted surfaces. Through innovative mechanical technology and human considerations, we have made C series to be quiet, safe, maintenance-free and cost-effective vacuum pumps.

Features

High chemical resistant

All wetted parts of C series pump are made of PTFE which are ideal for extremely aggressive / corrosive gases and vapors.

Long-term Durable

Head cover and diaphragm made of PTFE/PTFE-coated with stability core for unsurpassed long-term performance.

No air pollution, maintenance free

Chemvak C series pumps are driven by Diaphragm, without the need of lubricant, regular oil changes and maintenance; with no oil pollution.

Practical

> smooth surfaces for easy cleaning (NT)
> sealing system provides reduced leakage rates for improved ultimate vacuum

Quiet and low vibration

Driven direct by motor with no additional belt-driven transmission; the quality vibration-proof assembly makes Chemvak C series run at the lowest noise level among all other equivalent pumps.

Thermal protection device

Every motor of Chemvak C series pump has a built-in thermal protection device to shut off the pump automatically when overheated and then resume working when the temperature cools down.

International safety certification

CE certification

The stability core principle: for unprecedented long-term performance

chemistry diaphragm pumps provide optimum performance and unsurpassed service intervals even in harsh chemical applications. We achieve this unmatched reliability by manufacturing the most highly stressed components – the head cover and clamping disk.

> high quality carbon-fiber-reinforced fluoroplastics provide long term chemical resistance
> this thick-walled, diffusion resistant, molded fluoroplastic is supported by a stable metallic core for durability
> mechanical precision finishing ensures reproducible WIGGENS quality

Chemistry diaphragm pumps, chemistry vacuum systems and chemistry pumping units


Our chemistry diaphragm pumps are available in a full range of volume flow rates and ultimate vacuum options. Connecting pump heads in series as two-stage pumps improves an ultimate vacuum to as much as 0.6 mbar. Connecting heads in parallel provides higher flow rates. Our line offers combinations that satisfy virtually any laboratory need.

> + convenient, quick and simple to use due to the side-mounted power switch

> + whisper quiet and very low vibration
> + requires minimal bench top space
> + maintenance-free drive system and proven long diaphragm life
> + high chemical resistance

C series diaphragm pump offers a compact, high performance and easy-to-use solution which is perfect for single filtrations. This development, based on the highly successful technology of the oil-free diaphragm pump, provides a well-proven and extraordinarily long diaphragm life time. The functional, space saving and innovative design with easily accessible, side-mounted power switch ensures convenient and quick operation for day-to-day lab work. Robust PTFE coated diaphragms and valves provide optimal chemical resistance. A manual vacuum regulator valve for adjustment of pumping speed and ultimate vacuum is available as accessory.

Chemical Resistant Diaphragm Pumps Flow curve

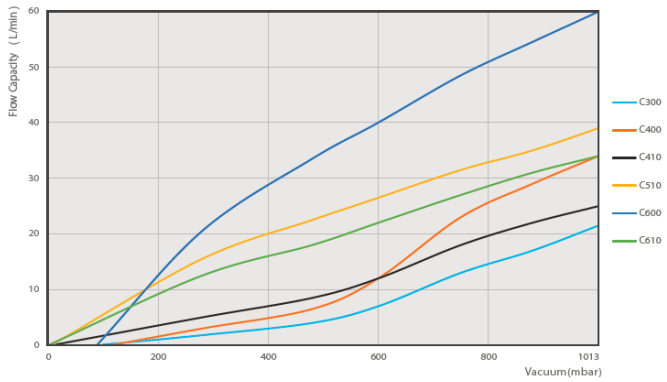Note: Above data is based on the 220V/50Hz instrument

Package 1

The pump with chemical resistant vacuum regulator (order no. 169311-06)

Package 2

The pump with VDM vacuum digital meter (order no. 169312-06)

Chemical Resistant Vacuum Controller DVR480

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WIGGENS China
Room 303, Hall C, Office Building M8, No.1 Jiuxianqiao East Road,
Chaoyang District, Beijing 100015, China
Tel: +86 400-809-2068
Fax: +86 400-809-2068-112
Email: info@ wiggens.com service@wiggens.com
Website: www.wiggens.com

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WIGGENS GmbH